

SENCOR®

SOC 4200BL / SOC 4201GD



TOOTHBRUSH WITH INDUCTIVE CHARGING



Translation of the original manual



SENCOR®

SOC 4200BL / SOC 4201GD



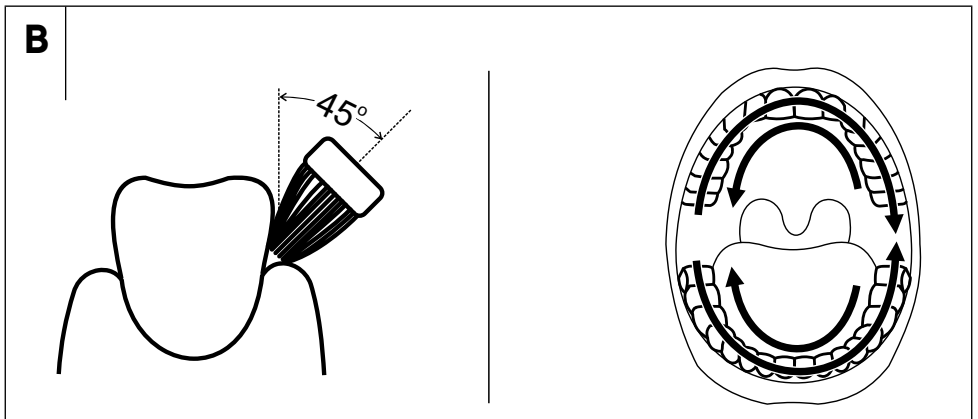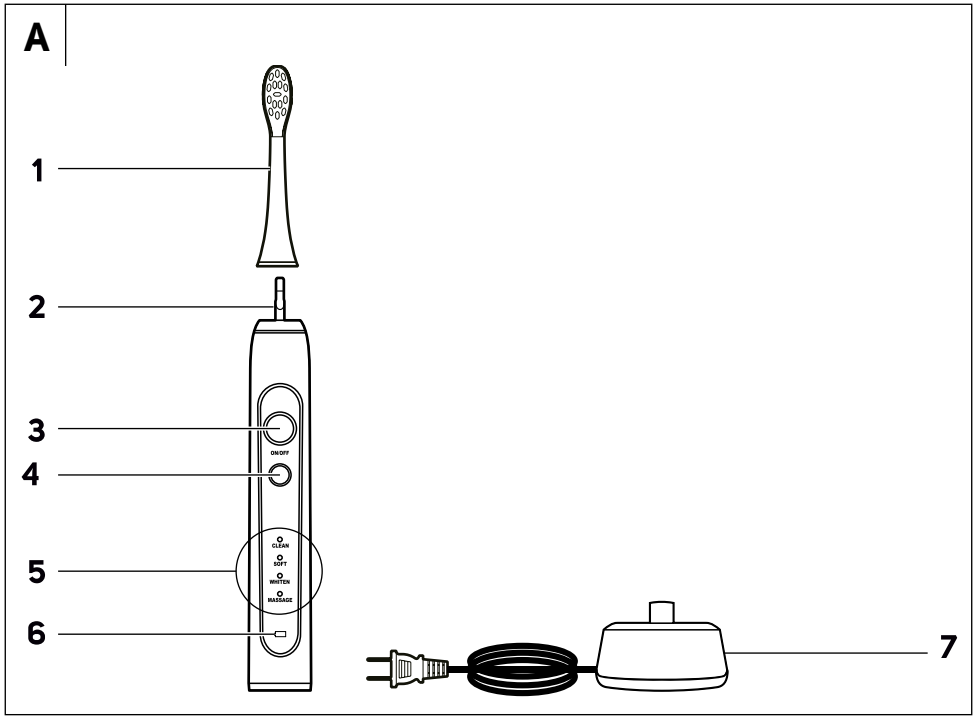

EN Toothbrush with inductive charging

Important safety instructions

READ CAREFULLY AND STORE FOR FUTURE USE.

GENERAL WARNINGS

- Toothbrushes may be used by children and persons with reduced physical or mental capabilities or with a lack of experience and knowledge, if they are supervised or have been instructed to use the appliance in a safe manner and understand the potential dangers.
- Children must not play with the appliance. Cleaning and maintenance performed by the user must not be performed by unsupervised children.
- In the event that the power cable is damaged, have it repaired at a professional service centre to prevent causing a dangerous situation. It is forbidden to use the appliance if it has a damaged power cable.
- The parts of the appliance that need to be mounted or stood upright, must be mounted or stood upright in such a way that they cannot fall into water.
- Use the supplied original charging base to charge the toothbrush. Before connecting the charging base to a power socket, make sure the voltage stated on its rating label corresponds to the voltage in your power socket.

ELECTRICAL SAFETY

- Before connecting the charging base to an el. power socket, make sure the voltage stated on its rating label corresponds to the voltage in your power socket.
- Connect the appliance only to a properly grounded socket. Do not use an extension cord.
- Do not connect and disconnect the power cord of the charging base to and from the power socket with wet hands.
- Ensure that the power plug contacts do not come into contact with water or any other liquid. If precisely this happens, carefully dry the power plug before using it again.
- Do not disconnect the charging base from the power socket by pulling on the power cord. This could damage the cord or the el. power socket. Disconnect the power cord from the power socket by pulling on the plug.
- Do not place heavy items on the power cord. Make sure the power cord does not hang over the edge of a table and that it is not touching a hot surface or sharp objects.
- Keep the power cord dry.
- If the power cord, the USB power cable or the USB adapter is damaged, it must be replaced by the manufacturer, its service technician or by a similarly qualified person, as this will prevent a dangerous situation from arising. It is forbidden to use a damaged power cord, USB power cable or USB adapter.
- Do not wind the power cord around the charging base.
- The appliance is equipped with a built-in long-life battery that can only be removed at a professional service centre. Never attempt to remove it from the appliance yourself.



Attention:

Never disassemble the appliance for the purpose of removing the battery. Never short circuit the appliance or throw it into a fire. Incorrect use, risks releasing toxic substances or causing a fire or explosion.

- Do not repair the appliance yourself or make any adjustments to it. Have all repairs or adjustments of the appliance performed at an authorised service centre. Tampering with the appliance during the warranty period may void the warranty.
- To prevent the risk of injury by electrical shock, never submerge the charging base and its power cord, the USB charging cable and USB adapter in water or any other liquid, or wash these parts under running water.
- If the charging base falls into water, under no circumstances continue using it, but rather hand it over to an authorised service centre for inspection.

SAFETY DURING USE

- This toothbrush serves to clean teeth and gums in the household. Do not use it for commercial purposes or for purposes other than those for which it is intended.
- Do not expose the toothbrush or its accessories to the effects of the outdoor environment, extreme temperatures, direct sunlight or locate it in the vicinity of an open flame, heating elements and other equipment or appliances that could be potential sources of heat.
- Always locate the charging base on an even and stable surface.
- Keep the power cord plug of the charging base dry.
- Charge, use and keep the toothbrush at room temperature.
- For hygiene reasons, toothbrush heads should be used by only a single user. Regularly replace the toothbrush head every 3 months, or sooner if it shows signs of wear. Use only original replacement toothbrush heads.
- When the toothbrush is used by children, supervise them while they are brushing their teeth until they acquire the proper tooth brushing technique.
- In the event that dental work or oral surgery was performed in the past two months, consult the use of a sonic toothbrush with your dentist.
- In the event that gum bleeding occurs when the toothbrush is used and if bleeding does not stop even after a week of its use, contact your dentist.
- If your toothpaste contains peroxide, sodium bicarbonate or another bicarbonate (frequently used in whitening toothpastes), then thoroughly clean the brush and the handle after use.
- In the event that using the toothbrush causes unpleasant or painful feelings, stop using the toothbrush and consult this situation with your dentist.
- In the event that you have a pacemaker or another implanted device, consult the use of the toothbrush with your doctor or manufacturer of the implant.
- Perform cleaning and maintenance according to the instructions in chapter Cleaning and maintenance. Do not wash the toothbrush or its accessories in a dishwasher.
- Do not make any modifications to the toothbrush or its accessories. By tampering with the product you risk voiding your legal rights arising from unsatisfactory performance or quality warranty.

- Thank you for purchasing a SENCOR brand product and we hope that you will be happy with it.
- Prior to using this appliance, please read the user's manual thoroughly, even in cases, when one has already familiarised themselves with previous use of similar types of appliances. Only use the appliance in the manner described in this user's manual. Keep this user's manual in a safe place where it can be easily retrieved for future use. In the event that you hand this appliance over to somebody else, make sure to also include this user's manual.
- Carefully unpack the appliance and take care not to throw away any part of the packaging before you find all its parts. We recommend saving the original cardboard box, packaging material, purchase receipt and responsibility statement of the vendor or warranty card for at least the duration of the legal liability for unsatisfactory performance or quality. In the event of transportation, we recommend that you pack the appliance in the original box from the manufacturer.

we recommend charging it at regular intervals (e.g. after 2-3 weeks of use).

DESCRIPTION OF THE APPLIANCE

- | | |
|--|-----------------------------|
| A1 Toothbrush head with a coloured ring | A6 Battery indicator |
| A2 Connecting pin | A7 Charging base |
| A3 Main switch | |
| A4 Mode toggle switch | |
| A5 LED mode indicators: | |
| CLEAN – standard cleaning mode | |
| SOFT – gentle cleaning mode for sensitive gums | |
| WHITEN – whitening mode | |
| MASSAGE – gum massaging mode | |



Note:

The toothbrush remembers the last used mode and therefore it is not necessary to reactivate it by pressing the mode toggle button.

If you press the mode toggle button while cleaning your teeth, the toothbrush will switch to the new mode and the respective light will be lit.

- Place the toothbrush bristles on to the teeth at an angle of 45° in line with the gum, see figure B. Using short slow motions, move the toothbrush to the other teeth. Once the outer side of the upper right dental arch is brushed, move the toothbrush to the inner side of the dental arch. Again, proceed from back to front. Adjust the toothbrush angle to the tooth arrangement. When brushing the inner side of the front teeth, hold the toothbrush handle vertical so that the tip of the toothbrush head points upwards. Finally, brush the upper tooth surfaces.
- In the same way, brush the right half of the lower dental arch, the left half of the lower dental arch and finally also the left half of the upper dental arch. When brushing the inner side of the front teeth of the lower dental arch, hold the toothbrush handle vertical so that the tip of the toothbrush head points downwards.

BEFORE FIRST USE

- Before first use, rinse the heads of the toothbrush with clean warm water.
- Wipe the body of the toothbrush with a dry or a very lightly dampened wiping cloth. Never submerge the charging base and its power cord in water, do not clean these under running water or wash them in a dishwasher.

CHARGING

- Place the charging base on an even and stable surface and connect the power cord to a power socket.
- Place the turned off toothbrush on the charging base. Charging will commence automatically and the battery LED indicator on the body of the toothbrush will light up.
- When the battery is fully charged, the LED indicator will turn off.



Note:

When charging the toothbrush for the first time after purchase or when it has not been used for more than three months, it is necessary to charge it for 24 hours. The toothbrush may heat up during the charging process. This is normal and does not indicate a malfunction.

- When the charging process is complete, we recommend disconnecting the charging base from the power socket.
- The toothbrush is designed to regularly run for up to 3 weeks after being charged.
- Although it is possible to recharge the toothbrush after it is fully discharged, in order to maintain the optimal lifetime of the battery,



Note:

The automatic shut-off timer automatically turns the toothbrush off after 2 minutes of operation. This is the tooth brushing time recommended by dentists for the sonic toothbrush.

- When you have finished using it, turn off the toothbrush using the main switch and clean it out according to the instructions contained in the chapter Cleaning and maintenance. Then place it into the storage glass or into the case.

CLEANING AND MAINTENANCE



Warning:

Prior to cleaning, the toothbrush must always be turned off and the charging base must be disconnected from the el. power source.

For cleaning the toothbrush, never use cleaning products with an abrasive effect, solvents etc. Otherwise, the surface of the appliance may be damaged.

- For hygiene reasons, regularly clean the toothbrush after every use.
- Thoroughly rinse the toothbrush head under lukewarm running water to remove any toothpaste or other remains. If whitening pastes are used when brushing teeth, it is recommended to apply a drop of soap on to the toothbrush head and then to thoroughly wash it under running lukewarm water. Remove the toothbrush head from the handle at least once per week and clean it under running water in the location where it is attached to the handle.
- Rinse the toothbrush handle under running lukewarm water and then dry it.
- In the event that the charging base is soiled, wipe it with a clean dampened cloth and then wipe it dry. Ensure that the power plug contacts do not come into contact with water or moisture.
- Never dry parts of the toothbrush over heat sources such as wood stoves, electric or gas stoves, etc.

REPLACING THE TOOTHBRUSH HEADS

- Regularly replace the toothbrush head every 3 months, or sooner if it shows signs of wear. Use only original replacement toothbrush heads. Replacement heads are available from you vendor both in online and brick&mortar stores.

STORAGE

- When not using the toothbrush for an extended period of time, disconnect the charging base from the power socket and put the cleaned toothbrush and toothbrush head in a dry, clean place. Do not expose the stored toothbrush to extreme temperatures, direct sunlight and excessive humidity and do not store it in an overly dusty environment.

TECHNICAL SPECIFICATIONS

Toothbrush


Battery	Li-Ion / 3.7 V / 650 mAh
Operating time when fully charged	up to 100 minutes
Noise level	65 dB(A)
Protection class level against damaging entry of water	IPX7

Charging base

Input	100-240 VAC, 50/60 Hz, 0,8-1,6 W
Electrical shock protection class	II
Protection class level against damaging entry of water	IPX7

The declared noise emission level of the appliance is 65 dB(A), which represents a level A of acoustic power with respect to a reference acoustic power of 1 pW.

Explanation of technical terminology

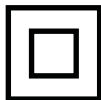

Safety class for electrical shock protection:
Class II – Electrical shock protection is provided by double or heavy-duty insulation.

Protection class level against damaging entry of water:
IPX7 – The appliance is protected against damaging entry of water when submerged at a depth of 1 metre for 30 minutes.

Changes to text and technical parameters are reserved.

INSTRUCTIONS AND INFORMATION REGARDING THE DISPOSAL OF USED PACKAGING MATERIALS

Dispose of used packaging material at a site designated for waste in your municipality.

DISPOSAL OF USED BATTERIES

Prior to handing over the appliance at a specialised collection facility, it is necessary to ensure that the built-in battery is removed. Battery removal will be performed at the authorised service centre. In the event that the product has been disassembled for the purpose of ecological disposal then it must not be reassembled.



This symbol indicates that the used battery does not belong in ordinary municipal waste. Hand over the used battery to designated collection points for proper disposal. Correct disposal of the battery will help prevent negative impacts on the environment and human health.

Disposal in other countries outside the European Union.

This symbol is valid in the European Union. If you wish to dispose of this product, request the necessary information about the correct disposal method from the local council or from your retailer.

INSTRUCTIONS AND INFORMATION REGARDING THE DISPOSAL OF USED PACKAGING MATERIALS

Dispose of used packaging material at a site designated for waste in your municipality.

DISPOSAL OF USED ELECTRICAL AND ELECTRONIC EQUIPMENT



This symbol on products or original documents means that used electric or electronic products must not be added to ordinary municipal waste. For proper disposal, renewal and recycling, hand these appliances over to determined collection points. Alternatively, in some European Union states or other European countries you may return your appliances to the local retailer when buying an equivalent new appliance. Correct disposal of this product helps save valuable natural resources and prevents potential negative effects on the environment and human health, which could result from improper waste disposal. Ask your local authorities or collection facility for more details. In accordance with national regulations penalties may be imposed for the incorrect disposal of this type of waste.

For business entities in European Union states

If you want to dispose of electrical or electronic equipment, ask your retailer or supplier for the necessary information.

Disposal in other countries outside the European Union.

This symbol is valid in the European Union. If you wish to dispose of this product, request the necessary information about the correct disposal method from the local council or from your retailer.



This product meets all the basic requirements of EU directives related to it.

